Smilyss®

Ready to Smile

by **SPORT-ELEC**
INSTITUT

FR

ES

NL

DE

GB

RO

Mode d'emploi

Jet dentaire

Smilyss

SPORT ELEC® Institut - 27520 BOURGTHEROULDE - FRANCE

www.sport-elec.com

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SAFETY INSTRUCTIONS

IMPORTANT SAFETY INSTRUCTIONS:

Before using your water jet please pay particular attention to the instruction manual, it contains important information on the product use, safety, care and maintenance.

Please store this instruction manual in a safe place so you can refer to the manual if you need to.

It is advisable to supply the device to operated circuit-breakers without integral and a rated power releasing no more than 30mA.

Only use this device with alternating current 230V. Before plugin, please check the voltage of your home matches the voltage device. You can refer to the device identification sheet.

The device must be plugged in a safe standards electrical distribution connection.

CAUTIONS:

NEVER IMMERSER THE APPLIANCE IN WATER OR ANY OTHER LIQUID.

DO NOT USE ABRASIVE CLEANERS OR HOUSEHOLD SOLVENTS.

GB

- If the appliance stops unexpectedly during the use, failed or if the device or its cord are damaged, switch off and unplug the device.
- Always unplug the device from the electrical outlet before cleaning or maintenance, and after use.
- Do not let appliance cords hang from counter-tops.
- Keep the appliance and cable away from sources of heat.
- Always use a suited and all working extension cord.
- Do not pull the cord to remove the appliance from the outlet.
- Keep the appliance away and out of reach of children.
- Do not use in wet environments (sauna, hammam) and near electrical appliance or trunk main with earthing connection.
- Never immerse the appliance in water or any other liquid.
- Do not use outdoors. This appliance has been designed for domestic indoor use only.
- Do not put the device down when it is in motion.
- Only use the accessories provided with the

appliance.

- Switch off the power supply and unplug from the wall before cleaning and or changing accessories.
- Do not operate your appliance in the presence of explosive or smoke.
- Any intervention or repair must be done by an authorized repairer: unauthorized repairer may represent a real danger for the user.
- This appliance has been designed for domestic indoor use only. Inadequate or improper use of the device will cancel the warranty and won't commit the liability of the manufacturer.

CONNECTING

The cord of the base can be plug in any alternating current 230 volts (~) including the shaver outlets. The power pack may become hot when charging. It is not affecting its functioning, overload is excluded.

MAINTENANCE

Unplug the device from the outlet. Rinse the sprinkler head and the handle separately during a few seconds.

Caution: make sure the water doesn't go in the handle.

Put the head away in the compartment after cleaning. Put the handle back on the base. Regularly clean the charger base with a soft cloth. The base has to be perfectly dry before plugging into the outlet.

INSTRUCTIONS MANUAL

1 - Take off the tank. Make sure you're holding the base (1) from one hand, and the tank with the other hand. Fill the tank with warm water and put it back on the device.

2 – Plug the device in the outlet.

Note: before the first use, rinse the device allowing the water to pass one full water tank. If the device is used with a "mouthwash", you need to rinse the device with pure water after each use in order to remove any leftover. If the device has not been used for a while, it is advised, for sanitary conditions, to follow the instructions for initial use and allow a full water

container to pass.

3. Choose a water jet head or a multi jet head and place it on the base handle, this way the colored head is in the same direction as "the water stop" button.

4. Put the pressure control button at the minimum (in front of LO) for the first use. You can gradually increase the pressure in


order to find the most comfortable position (MI then HI/ turbo).

5. Put the device in motion by pushing the "on/off" button located on the base, lean over your sink and put the head brush in your mouth. Put the button "on/off", on the handle of the water jet, on the "on" (I) position.

Note: keep your mouth lightly open and your elbow raised so that the water can flow into the sink.

6. Run the water along your gums and between your teeth. Never point la sprinkler head directly at your gums. For best results, it is advisable to start with your back teeth and work forward the front teeth.

7. To pass from the "jet" function to the hydro massage spray (or vice-versa), use the multi jet head brush and run the colored part. Use the jet

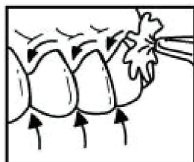for an accurate/ precise cleaning and a larger jet for a massage and cleaning in general.

8. Once you're done, put the "on/off" water jet (on the handle) on the "off" position (O). And turn off the unit by pressing the button on / off (on the base unit). Take away the head and put it away in the compartment. Empty the device from all liquid. Unplug the outlet from the wall.

Note: never let the device on the "on" position when

not in use. The water pressure will keep on form and may damage the device. Always set the "on/off" button from your water jet and the base on the "off" position (O) when the device is not in use.

You can also use the interdental brush head with integrated water jet) which combines brushing and jet efficiency. This tool allows a good oral hygiene, even in the hard to reach area.

ENVIRONMENT FRIENDLY DISPOSAL

Your appliance contains valuable materials which can be recovered or recycled. Please leave it at a collection point as listed by your local authority.

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