## GARANTIE

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## WARRANTY

This appliance is guaranteed for 12 months from the date of purchase. During the period of the guarantee, the distributor will take responsibility for defects occur, the appliance will be repaired or exchanged as the distributor deems necessary. Damage arising from improper use or normal wear and defects, which do not materially affect the performance of the appliance, is not covered under this guarantee. The guarantee is not valid if unauthorized repairs have been attempted or if replacement parts used are not approved by the supplier. The guarantee is valid only if it is stamped and signed by the retailer at the time of the purchase. For all claims made during the period of the guarantee, send the appliance and the warranty card to your retailer. All other claims are excluded from the guarantee unless specifically prohibited by law.

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Adresse /Anschri /  
Adress / Direcci n / Adres .....

Code Postal /Postleitzahl /  
Ville .....

Postcode / C digo Postal / Postcode  
Place .....

CACHET DU REVENDEUR / STEMPEL VOM H NDLER /  
STAMP OF RETAILER / SELLO DEL VENDEDOR /  
STEMPEL VAN DE VERKOPER

DATE D'ACHAT / EINKAUFSDATUM / DATE OF PURCHASE/  
FECHA DE COMPRA / VERKOOPSDATUM

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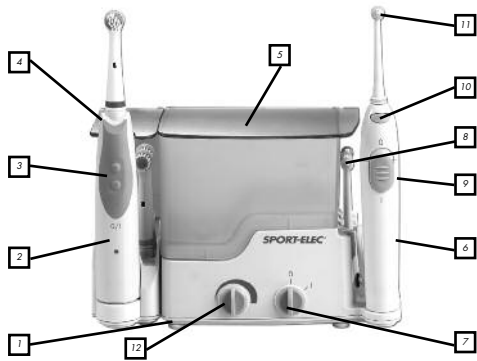
- 1 - Base
- 2 - Toothbrush
- 3 - Toothbrush handle with ON/OFF button
- 4 - Storage for 4 brush heads
- 5 - Water beaker
- 6 - Oral irrigator
- 7 - ON/OFF button of the irrigator's handle
- 8 - Storage for 4 nozzles
- 9 - ON/OFF button of the water jet
- 10 - Heads locking button
- 11 - Water jet control : direct stream and spray hydro-massage
- 12 - Water jet pressure control

#### NOTE OF USE

#### PRECAUTIONS BY USING AN ELECTRIC DEVICE :

When using electrical appliances, basic safety precautions should be always followed to avoid any accident. Read attentively all instructions before using the appliance, it's cleaning, charging and maintenance.

- Before starting your device read carefully the instructions manual. It contains important information concerning the use, the safety and the maintenance of your apparatus.

- It is recommended to use this appliance on an electrical installation using a differential with nominal current loss not exceeding 30 mA.

- While reading this manual, please, refer to explanatory pictures.

- Use this device only with AC current 230V. Before connecting the device make sure that the tension of your electrical network corresponds to the one mentioned on the bottom of the charge base. The device should be connected only to that electrical supply network installation of which is in conformity with the standards of safety.

- Switch off and disconnect your device in the following cases :

- in case of any anomaly during use, before any cleaning or maintenance, immediately after use.

- Do not let the cord hang on the edge of the table.

- Keep the cord away from heated and sharp surfaces

- Use an extension for your cord only if it is suitable and in a perfect state.

- Do not excessively bend or pull the cord. Unplug by grasping the plug, not the cord.

- Children do not recognize that electrical appliances can cause danger. Do not leave apparatus under operation within the range of children.

- Never submerge the apparatus in water or any other liquid.

- Do not place or store the device where it can fall into water ( tub, sink or toilet).

- Do not use the apparatus outdoors.

- Do not put your device anywhere when it is functioning.

- Do not use accessories other than those provided by the manufacturer with this unit to avoid danger. Unplug the device when you do not use it, before mounting or dismounting the accessories or before any cleaning.

- Clean the apparatus with a dry or slightly damp cloth. Don't use metal sponges or abrasive agents for cleaning.

- Make sure that the charge base is perfectly dry before plugging your device.

- Do not make the device function in presence of explosives or fumes.

- **WARNING :** To avoid any risk of fire or electric shock, do not dismantle the device - Any interference and repairing must be carried out by an approved specialist.

- Do not use your apparatus or unplug it or do not plug it if :  
- It is damaged.

- The electrical cord or the socket are damaged.

- You have noticed the abnormal functioning after some incident.

- Non-professional interferences can cause danger to the user.


- If the electrical cable is damaged, it should be replaced by the Manufacturer, its After-Sales service or a person with a similar qualification, in order to avoid any danger.

- This device is for household use only. Any inadequate use and non conformity to the instructions cancels the guarantee and cannot engage the responsibility of the manufacturer.

Don't use the apparatus in a damp room.

#### ANY DETEIORATION RESULTING FROM THE NON-OBSERVANCE OF THESE PRECAUTIONS OF USE SUPPRESSES AUTOMATICALLY THE GUARANTEE

#### SAVE THESE INSTRUCTIONS

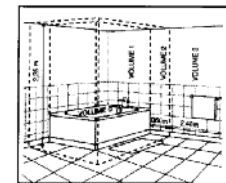This apparatus can only be installed in the bathroom in volume 3 (see the picture). It is prohibited to use it in volumes 1 et 2. It must be fixed at the wall using the screws, ankles and reception plate which come with it and not above the wash-hand basin.

#### CLEANING AND MAINTENANCE

Always empty the water beaker completely after use.

Note: if the unit is used with a mouthwash solution, a full water beaker should be run through by water after each use to rinse it and remove any traces of mouthwash solution.

#### TOOTHBRUSH :

To maintain the brushing head in good operating condition, carefully rinse it during several seconds under water. To rinse the toothbrush switch it on, it should operate. For a perfect hygiene we recommend soaking the brushing head in a denture cleaner once a week.

Clean the body of the toothbrush with a soft and slightly wet cloth. The toothbrush must be perfectly dry before placing it on its charging base.

#### Charging Base :

Dry the external sides of the charging base with a damp soft cloth.

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Note: Rinse the nozzles and the dental jet separately during a few. Attention: make sure that the water does not enter the handle.

3. Put back the nozzle into the storage after cleaning. Put back the handle on the base.

4. Clean the base with a soft and slightly wet cloth regularly. The base must be completely dry before plugging in the power socket.

#### ATTENTION :

**NEVER SOAK THE CHARGING BASE OF THE POWER SUPPLY UNIT IN WATER.**

**NEVER USE ABRASIVE PRODUCTS OR SPONGES**

#### MEDICAL SAFETY MEASURES FOR THE DENTAL JET

Please mind the following in order to guarantee your safety.

1. Consult your dentist before using the dental jet if you underwent an oral intervention or periodontal.

2. Consult your dentist in the case of an excessive bleeding of the gums. Do the same if this bleeding persists after 2 weeks.

#### HOW TO USE YOUR ORAL IRRIGATOR

1. Remove the beaker. Pay attention to hold well the base ( 1 ) in its place by your hand, whereas you raise the beaker by the other hand. Fill the beaker with warm water and put it back on the device.

2. Plug the device to outlet

Note : Prior to initial use, rinse the unit by letting water pass through a full water beaker. If the unit is used with a mouthwash solution, it should be rinsed with clean water after each use to remove any traces of mouthwash. If the unit has not been used for an extended period of time, to maintain hygiene it is recommended that you follow the instructions for initial use: let a full water beaker pass through.

3. Choose a nozzle (8) (there are 4 different colours for several users) and attach it on the handle so that its coloured head faces the same direction as water 'stop' switch (Picture 1)

4. Put the water jet pressure control on the base of the device into a minimum position (turn clockwise during the first use). You can increase the pressure gradually to find the most comfortable position.

5. Before plugging the device, make sure that ON/OFF button of the water jet of the handle is in the position 'off' (O). Start your dental jet by pushing the switch from the position "Off" to position « I » of the base, lean over the sink and put the nozzle in your mouth. Push the button of the water jet into position 'on' ( I ) .

Note : Keep your mouth slightly open and hold your elbow at height of your mouth, so that water can run in the sink.

6. Direct the water jet along your gum line and also between your teeth.

The nozzle should never be in direct contact with your gum.

For better results, it is advised to begin with molars towards the canines.

Consult your dentist for further information concerning use.

7. To change the function "direct jet" to the function "spray hydro-massage" (or vice versa), turn the coloured nozzle (Picture 9). Use the jet for a precise cleaning and the spray hydro-massage for usual cleaning and massage.

Note : In order to avoid splashes start the dental jet when it is out of your mouth. Hold the apparatus above the sink and point it downward.

8. Once finished, put the control switch to the position 'off' (O) and unplug the device. Remove the nozzle by pressing on the lock button and put it back into the storage. Empty the liquid of the apparatus. Unplug the outlet

Never leave the device in the position 'on' (I) when it is not used. The pressure of water will continue to form and can damage the device. Always put ON/OFF button of the water jet and the ON/OFF button in the position ' off ' (O) when the device is not used.

#### HOW TO USE THE TOOTHBRUSH

CONNECTION OF THE APPARATUS - CHARGING THE BATTERY  
The charging base can be connected to any plug of AC current 220-240 (-), including shaving sockets.

Place the toothbrush on the charging base. The green charging light will light up indicating the handle under charging. The base may warm up slightly during charge. This has no effect on the unit's operation. Overcharging is impossible.

Before the first use and in order to prolong the battery life, it is necessary to proceed as follows:

1. charge the battery during 16 hours ( full charge duration).

2. then discharge the battery completely by starting the toothbrush until automatic stop.

3. recharge to the maximum the battery of the toothbrush during 16 hours. This operation "discharging / recharging" is to be repeated every six months.

N.B. : When the apparatus starts to discharge, the intensity of brushing becomes definitely lower and automatically stops functioning. So it is necessary to place back the toothbrush on its base and to recharge the batteries.

After having carried out the operation "discharging / recharging" described above, you can use the toothbrush.

**NB :** To obtain maximum power of brushing we recommend keeping your brush on the charging base when not in use to maintain a full battery charge.

- Fix a toothbrush head to the supply block.  
- Before the first use of a toothbrush head, rinse it carefully under the water.  
- Apply your ordinary toothpaste.

- Rinse your mouth

- To avoid the splatters, apply the brushing head against your teeth before starting the unit.

- Place the brushing head on the surface of your teeth, then switch on the toothbrush and make circular motions on your teeth from top to bottom for the teeth best cleaning and gum massaging effect.  
It is not necessary to press strongly on the brush.

Clean all areas of your teeth including the inner and chewing surfaces and behind the back teeth.

- To eliminate the dental plaque completely, brush your teeth during two minutes at least.

- Rinse your mouth carefully after the brushing.

- When you first use the toothbrush, you may experience slight bleeding. This should stop after a few days when your gums get used to the toothbrush. If the bleeding persists during fifteen days we recommend you to consult your dentist.

The most effective method brushing consists in making gently pass the toothbrush from one tooth to the other, without pressing it.

The form of the brushing head change gradually within the time of use, what is normal and this also depends on your personal way of brushing your teeth. This deformation has no incidence on the efficiency of the brush. Do not however forget to change the head of brushing as regularly as a manual toothbrush.

- After use your toothbrush, rinse its brush head under tap water.

- Put your tooth brush back to its charging base.

#### WALL MOUNTING

Use only the mural support to ensure the maintenance of the apparatus and to avoid the falls into the sink.

Fix the screws and the dowels to the wall. Now you can shift the mouth shower out upward. Make the cord pass between the support and the charging base. Do not use the cord longer than necessary to reach the socket outlet.

#### Technical features :

Power supply: 230 V ~ 50 Hz - Intensity: 15 W

Batteries : 2,4 Volts. Operation life full charge 30-45 minutes.

Protection index : IP24 (index of protection against the solid foreign bodies larger than 12,5 mm and against the splashes and water projections.  
CE (89/336/CEE) EMC

#### ENVIRONMENTAL PROTECTION

Remove the rechargeable battery at the end of the unit's life.

Before removing, the battery must be completely emptied. Let the apparatus function until it stops alone.

Unscrew the bottom only to remove the battery. The apparatus will not be waterproof any more after having unscrewed the bottom.


La brosse à dents contient des batteries au Nickel Cadmium. The toothbrush contains Nickel Cadmium batteries. In a concern of environmental protection, you are asked not to throw your toothbrush to household refuse when it arrives at the end of life without having removed the battery from it beforehand.

When you throw your product at the end of its life, respect your national or local regulations for environmental protection and recycling of raw materials. We thank you for depositing it in a place provided for that purpose to secure its elimination in an environmentally respectful way.

#### RoHS

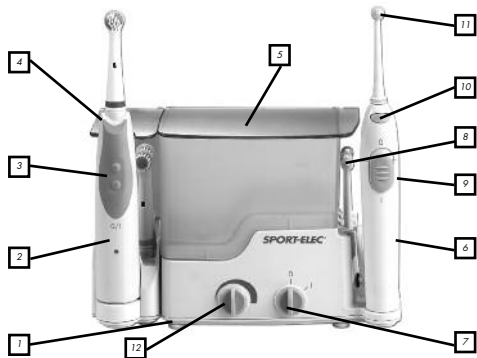(2002/95/EC) : This product is manufactured without materials which could have a harmful effect on ecology.

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